

B. Braun Vascular Systems

Uni-Graft® K DV

...the knitted impregnated vascular graft




Established product quality with a safety record of over 20 years

Uni-Graft® K DV


The knitted double-velour architecture of Uni-Graft® K DV, permit the inward migration of cells and ensures a firm integration into the surrounding tissue. The velours on the interior surface facilitates adhesion of the fibrin film which develops during the first phase of revascularisation.

The gelatine impregnation effectively seals all pores in the graft and therefore does not require preclotting. Furthermore, the impregnation prevents the penetration of blood into the graft wall.



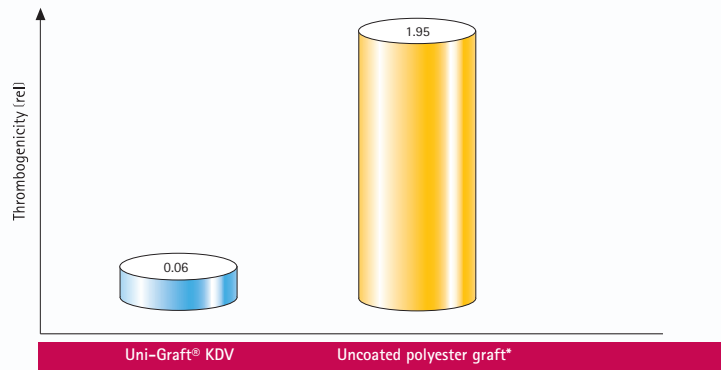

Uncoated polyester grafts (upg*) are by two orders of magnitudes more thrombogenic than Uni-Graft K DV

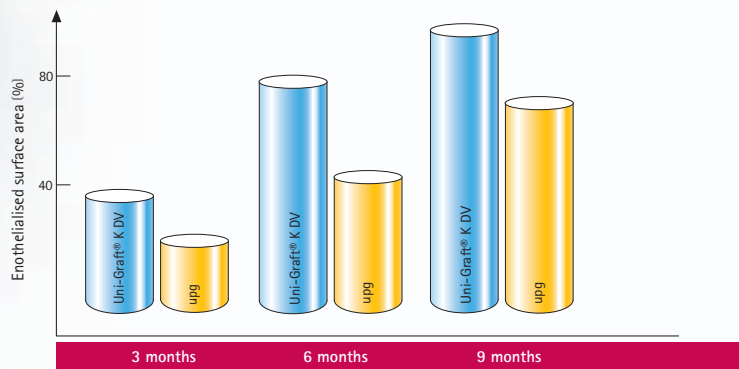

For these reasons the neointima develops more rapidly in the Uni-Graft® K DV graft as compared to in uncoated reference grafts. This was examined in a morphological study based on a canine model.

Ordering information

Uni-Graft® K DV – Bifurcated Prostheses

usable length	diameter	code number	Sales unit-pcs.
40 cm	12 x 6 mm	1104535	1
	12 x 7 mm	1104519	1
	14 x 7 mm	1104527	1
	14 x 8 mm	1104551	1
	16 x 8 mm	1104560	1
	16 x 9 mm	1104578	1
	18 x 9 mm	1104586	1
	18 x 10 mm	1104594	1
	20 x 10 mm	1104608	1
	20 x 11 mm	1104616	1
	22 x 11 mm	1104624	1
	24 x 12 mm	1104640	1

Uni-Graft® K DV – Straight Tubes

usable length	diameter	code number	Sales unit-pcs.
10 cm	6 mm	1101030	1
	10 mm	1101803	1
15 cm	6 mm	1101021	1
	7 mm	1101048	1
	8 mm	1101013	1
	10 mm	1101102	1
	12 mm	1101129	1
	14 mm	1101145	1
	16 mm	1101161	1
	18 mm	1101188	1
	20 mm	1101200	1
	22 mm	1101226	1
	24 mm	1101242	1
	26 mm	1102265	1
	28 mm	1102281	1
30 mm	1102303	1	
32 mm	1102320	1	
34 mm	1102346	1	
30 cm	6 mm	1102060	1
	7 mm	1102079	1
	8 mm	1102087	1
	10 mm	1102109	1
	12 mm	1101137	1
	14 mm	1104144	1
	16 mm	1104160	1
18 mm	1104187	1	

usable length	diameter	code number	Sales unit-pcs.	
	20 mm	1104209	1	
	22 mm	1104225	1	
	24 mm	1104241	1	
	26 mm	1102192	1	
	28 mm	1102206	1	
	30 mm	1102214	1	
	32 mm	1102222	1	
	34 mm	1102230	1	
	36 mm	1102249	1	
	40 cm	6 mm	1101269	1
		7 mm	1101277	1
8 mm		1101285	1	
10 mm		1101307	1	
18 mm		1101684	1	
60 cm	6 mm	1104063	1	
	7 mm	1104071	1	
	8 mm	1104080	1	
	10 mm	1104101	1	
	12 mm	1104128	1	
	70 cm	6 mm	1104421	1
7 mm		1104422	1	
8 mm		1104423	1	
100 cm	6 mm	1104403	1	
	8 mm	1104411	1	
	10 mm	1104420	1	

Uni-Graft® K DV – Tapered

usable length	diameter	code number	Sales unit-pcs.
60 cm	6 – 8 mm	1102869	1

Uni-Graft® K DV – Axillo-Bifemoral

usable length	diameter	code number	Sales unit-pcs.
100/60 cm	8 mm, right side-branch	1104705	1
	8 mm, left side-branch	1104713	1
	10 mm, left side-branch	1104730	1

Innovations for Vascular Surgery

...offers you what you expect from a safe and reliable vascular graft!



Indications:

- reconstructive procedures of the abdominal aorta and of peripheral vessels
- patients requiring systemic heparinisation
- emergency surgery
- patients with clotting disorders
- extra-anatomical procedures

Knitted double-velours:

- ensures good graft integration and tissue incorporation
- simple handling
- no fraying
- excellent biocompatibility

Gelatine impregnation cross-linked with isocyanate

- absolutely impermeable to blood
- no preclotting
- for immediate implantation
- athrombogenic
- accelerates cellular integration
- aldehyde free cross-linking

B. Braun Vascular Systems offers you competent services and high quality products for interventional procedures and vascular surgery:

- ▶ Silver Graft
- ▶ VascuGraft PTFE/SOFT
- ▶ Uni-Graft®
- ▶ ProVena
- ▶ Patches
- ▶ Insitucat
- ▶ Venostrip
- ▶ VascuFlex SE/SEC



For product information
please contact:

B. Braun Melsungen AG
Division Aesculap - Vascular Systems

Sieversufer 8
12359 Berlin
Germany

Phone +49 30 689897-0
Fax +49 30 689897-30

www.bbraun.com

B | BRAUN
SHARING EXPERTISE

Aesculap AG & Co. KG

Am Aesculap-Platz
78532 Tuttlingen
Germany

Phone +49 7461 95-0
Fax +49 7461 95-2600

www.aesculap.de

All rights reserved. Technical alterations are possible. This leaflet may be used for no other purposes than offering, buying and selling of our products. No part may be copied or reproduced in any form. In the case of misuse we retain the rights to recall our catalogues and pricelists and to take legal actions.